R5

Neck Care Massager

ELEEEELS

Instruction Manual

Clamping and
Kneading
Massage
Mechanism

Hot
Compress


Cordless and
Portable

2-Hour
Battery
Life


R5 at a Glance

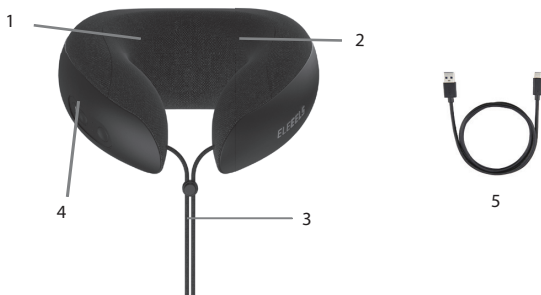- 1) Clamping and Kneading Massage
- 2) Hot Compress
- 3) Adjustable String
- 4) USB-C Charging Port
- 5) USB-C Cable

Specifications

- Rated Voltage: 7.4V
- Input: DC 5V/2A
- Battery Capacity: 2000 mAh
- Rated Power: 10W
- Weight: 1202 grams
- Operating Time: up to 1-1.5 hours of use per charge

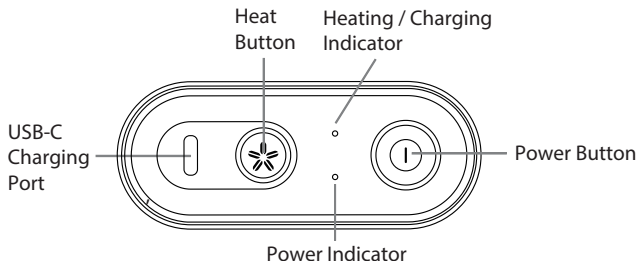INSTRUCTIONS

PLEASE READ PRIOR TO USE

CHARGING

- 1) Connect USB-C cable to the charging port.
- 2) The charging speed depends on what adapter is used. It is suggested to use an USB adapter supporting 5V/2A to fully charge in approximately 3 hours.
- 3) During charging, the charging indicator is in red solid light. When the charging indicator turns off, and the power indicator is in green solid light, the device is fully charged.

OPERATING

- **Put on the device. Adjust the string to the most comfortable position. Alternatively, lie on the sofa or chair, with the device supported by sofa or chair back. Relax your neck by leaning on the device.**
- Press the power button for 2 seconds to turn on the device. The default setting is Soothing Mode.
- Press the power button repeatedly to switch between Active Mode, Combination Mode, and Soothing Mode.
- Press the heat button once to turn on hot compress. Press the heat button once again to turn off hot compress.
- When the device stops operation after a 10-minute session, press the power button for 2 seconds to start a new one.
- Press the power button for 2 seconds to turn off the device.
- When the battery level is lower than 10%, charging indicator flashes for 5 times. Operation stops. Please recharge the battery.

• MAINTENANCE, CLEANING & STORAGE

- Use a slightly damp towel to clean the device and a soft cloth to dry it.
- Do not clean the device with chemicals or alcohol.
- Only use the USB-C cable supplied by the ELEELS for recharging.

WARNINGS & CAUTION

- Do not immerse in water. Keep away from liquids or heat sources.
- Do not remove screws or attempt to disassemble.
- Do not operate continuously for more than one hour. Allow device to rest for 10 minutes before re-using.
- For adult use only. Do not use if injured. Consult your doctor before use.

SAFETY INSTRUCTIONS

DANGER: TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

- Do not drop the device.
- Do not tamper or alter the device in any way.
- Never leave the device unattended during operating or charging.

DISCLAIMER

Please do not use the ELEEEELS R5, or any massage device without first obtaining approval from your doctor if you have any of the following conditions: severe cervical hyperosteoegeny, abnormal cervical vertebra, severe blood circulation disorder in the neck, neuropathy, pacemaker user, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, spondylosis, recent joint replacements, surgical pins or plates, or you have any concerns about your physical health.

WARRANTY

IF WITHIN ONE YEAR FROM THE DATE OF PURCHASE, THIS PRODUCT FAILS DUE TO A DEFECT IN THE MATERIAL OR WORKMANSHIP, ELEEEELS WILL REPAIR OR REPLACE THE PRODUCT; OR FAILED COMPONENT(S), FREE OF CHARGE

THIS WARRANTY EXCLUDES:

- (a) Damage caused by accident, abuse, or mishandling during transport;
- (b) Damage arising from use of units from unauthorized repair;
- (c) Malfunction of units that are not used in accordance with ELEEEELS instructions;
- (d) Damage for which the cost of repair exceeds the cost of the product;
- (e) Deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises, and
- (f) Any product for which a dated proof of purchase cannot be provided.

Please visit below link for instruction manual in other languages:
<http://www.eleeeels.com/im>

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