S2 | Hot Stone Massage Wand Collection


Instruction Manual

ELEEELS

Pro Massage Made Easy

ELEEELS S2 / S2A / S2B at a Glance

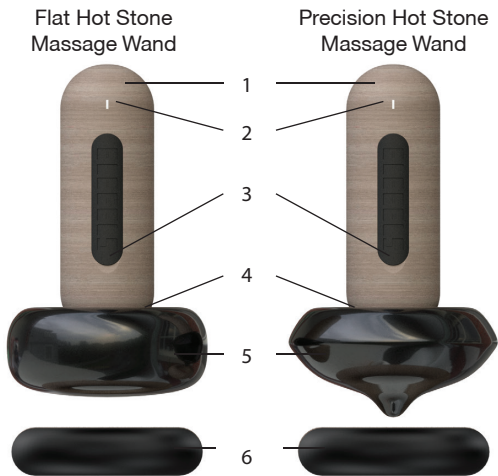7

8

- 1) The wands
- 2) Switch for temperature control
- 3) Type C charging port beneath the flip
- 4) Temperature indicator and power indicator
- 5) Hand polished natural bian massage stone
- 6) Torus holder
- 7) Aromatherapy massage oil (included in S2 only)
- 8) Type-C cable

Weight: 3.31 lbs. / 1.5 kg
Input: Type-C 5V/2A

CHARGING


- 1) Flip over the rubber cover to unveil the USB Type-C charging port. Connect USB Type-C cable to the charging port.
- 2) The charging speed depends on what adapter is used. It is suggested to use an adapter supporting Type C (5V/2A) to fully charge in approximately 4 hours.
- 3) During charging, the power indicator flashes in green slowly. When the green light stops flashing, the device is fully charged.
- 4) The power indicator shows the battery level:
 - a) Flashing Red Light: 0-10%
 - b) Flashing Green Light: 10-99%
 - c) Solid Green Light: 100%

OPERATING - HOT STONE MASSAGE WAND


- 1) Switch the control clockwise once to power on the device. Switch the control repeatedly to adjust the temperature setting during operation. Switch clockwise for higher temperatures, and anti-clockwise for lower temperatures or power off.
- 2) The temperature indicator indicates the corresponding heat level.
 - Level 1 (yellow light, 45 degrees Celsius)
 - Level 2 (orange light, 50 degrees Celsius)
 - Level 3 (red light, 55 degrees Celsius)
 - Level 4 (white light, 60 degrees Celsius)
- 3) Always start with level 1 for first-time use.
- 4) After switching on for 1 minute, the wand(s) will be warmed up and ready to use.
- 5) Place the hot stone wands either on arms, shoulder, neck, back, abdominal, buttock, legs, or feet of your partner or yourself and apply the 7 massage skills to enjoy hot stone massage treatment.
- 6) Apply a small amount of massage oil on the skin or on the wand when lubrication is needed.
- 7) Switch the control anti-clockwise once during Level 1 to power off the device.
- 8) When the battery level is lower than 10%, the power indicator flashes in red. The device will turn off in 3 minutes. Please recharge the battery.

OPERATING - 7 MASSAGE SKILLS BY THE 2 WANDS

Flat Hot Stone
Massage Wand

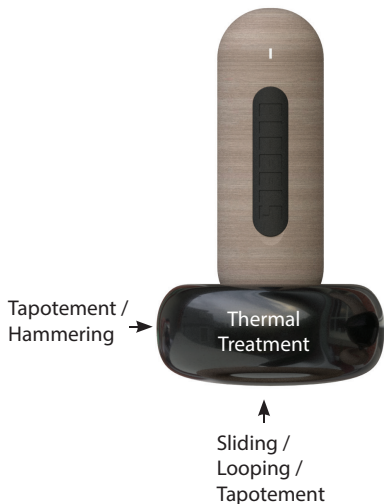Precision Hot Stone
Massage Wand

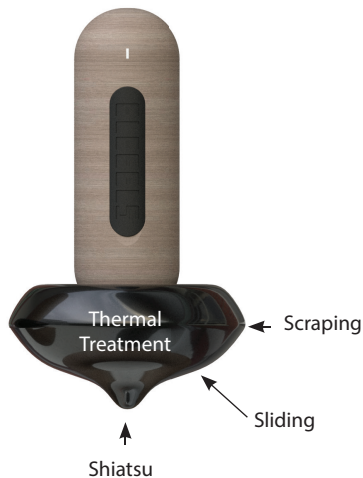Please visit below link for tutorial of massage skills by ELEEELS S2 / S2A / S2B:
<https://www.youtube.com/c/ELEEELS>

WARNINGS & CAUTIONS

- This device is specifically intended for body massage using hot bian stone(s).
- This device and its components are not designed to be used by children.
- Gently place the wand(s) on the torous holder.
- **Always touch the stones first to feel the temperature before applying on body.**
- Stop the massage immediately if you feel uncomfortable.
- Do not use the device inside Finnish bath with very high humidity.
- Always switch off the product when not in use and before cleaning.
- Avoid touching the charging port or power switch with wet or damp hands.
- Do not place the device next to a heat source, water, or other liquids.
- Do not place any objects or pour liquid on the device.
- If this product stops working, switch off the product immediately. Seek professional advice before attempting to re-use it.
- Do not remove screws or attempt to disassemble.
- Do not drop the appliance.
- Do not tamper or alter the appliance in any way.

MAINTENANCE & CLEANING

- Use a soft towel to clean the device. Add a minimal amount of detergent to clean the massage oil.
- Do not clean the device with chemicals or alcohol.

DISCLAIMER

Please do not use the ELEEELS S2 / S2A / S2B without first obtaining approval from your doctor if you are pregnant, or diagnosed with any of the following conditions: malignant tumor, scoliosis or abnormality, hypertension, osteoporosis, thrombosis, severe aneurysm, dermatitis and other skin infections (including inflammation of subcutaneous tissue), wounds, and severe circulatory disorders.

Our hot stone is a natural product and has been hand polished as such contains slight variations in color, veining, texture, contour, or shading.

WARRANTY

IF WITHIN ONE YEAR FROM THE DATE OF PURCHASE, THIS PRODUCT FAILS DUE TO A DEFECT IN THE MATERIAL OR WORKMANSHIP, ELEEEELS WILL REPAIR OR REPLACE THE PRODUCT; OR FAILED COMPONENT(S), FREE OF CHARGE.

THIS WARRANTY EXCLUDES:

- (a) Damage caused by accident, abuse, or mishandling during transport;
- (b) Damage arising from use of units from unauthorized repair;
- (c) Malfunction of units that are not used in accordance with ELEEEELS care instructions;
- (d) Damage for which the cost of repair exceeds the cost of the product;
- (e) Deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises, and
- (f) Any product for which a dated proof of purchase cannot be provided.

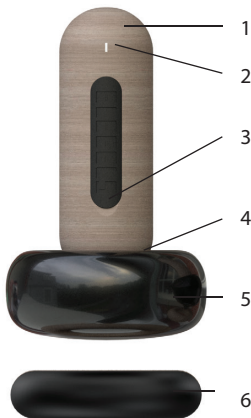
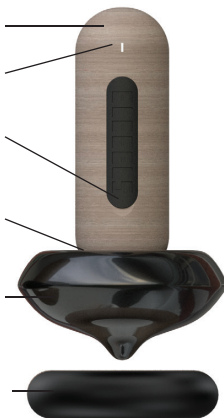
Please visit below link for instruction manual in other languages:

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