X1T

Percussive
Massage
Device

Instruction Manual

CORDLESS / RECHARGEABLE


24V Lithium-ion Battery

LIGHTWEIGHT DESIGN

2.4 lbs. / 1.1 kg


OPERATING

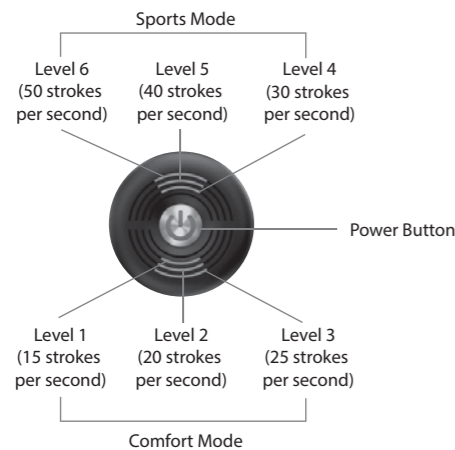- Press power button for 2 seconds to switch on.
- Press power button again for comfort mode level 1.
- During operation, press power button repeatedly for higher speed level. The highest speed level is sports mode level 6. The blue light speed level indicator indicates the corresponding speed level.
- During sports mode level 6, press power button to pause the device.
- Massage the desired body part at your preferred speed by applying pressure as needed.
- When the battery level is lower than 10%, the power indicator flashes in red. Operation stops. Please recharge the battery.
- When the massage device stops operating after a 10-minute session, press the power button to start a new session.
- When the device stops operation for 2 minutes, the device shuts off automatically. Alternatively, press power button for 2 seconds to switch off.

We developed the ELEEELS percussive massage device to accelerate recovery and prevent muscle damage in order to relieve pain and promote circulation.

ELEEELS X1T

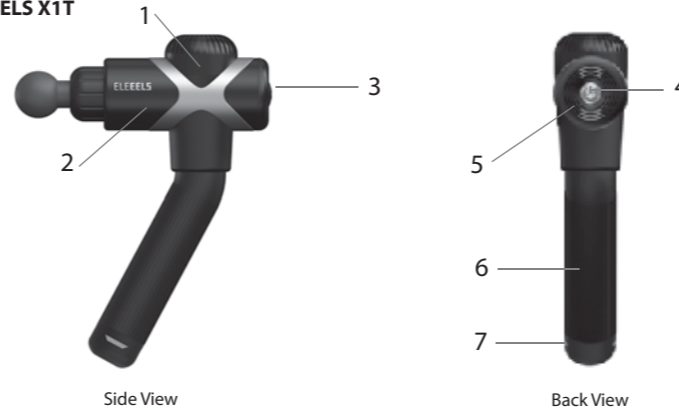- 1) High-torque, high-speed brushless motor (delivers up to 50 strokes per second)
- 2) Endurance Motion Technology™ (efficient energy use and long operation time)
- 3) Quiet Motion Technology™
- 4) Sports Mode Indicator
- 5) Comfort Mode Indicator
- 6) Rechargeable lithium-ion battery (up to 10 hours of use per charge)
- 7) Cordless and Lightweight (1.1 kg)

Vibration Frequency:
Comfort Mode: Level 1 (15 strokes per second), Level 2 (20 strokes per second), Level 3 (25 strokes per second)
Sports Mode: Level 4 (30 strokes per second), Level 5 (40 strokes per second), Level 6 (50 strokes per second)

Weight: 2.4 lbs. / 1.1 kg
Charging Port: USB Type-C

HOW DO I ADJUST THE ANGLE OF THE HANDLE?

Twist the battery (handle) 180 degrees until you hear a 'click'.

MAINTENANCE & CLEANING

- Use a slightly damp towel to clean the device and a soft cloth to dry it.
- Do not clean the device with chemicals or alcohol.
- Only use the USB Type-C cable supplied by the ELEEELS for recharging.

WARNINGS & CAUTIONS


- Do not immerse in water. Keep away from liquids or heat sources. Keep ventilation ports free from dust and debris.
- Do not remove screws or attempt to disassemble.
- Do not operate continuously for more than one hour. Allow device to rest for 10 minutes before re-using.
- For adult use only. Do not use if injured. Consult your doctor before using this product.

SAFETY INSTRUCTIONS

DANGER: TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

- Apply to a clean, dry clothed area of the body by lightly pressing and moving across the body for approximately 60 seconds per area.
- Only use the device on the body's soft tissue as desired. Do not use on head or any hard, bony area of the body. Discontinue use immediately if experiencing any pain.
- Do not drop the device.
- Do not tamper or alter the device in any way.
- Never leave the device unattended during operating or charging.

4 Interchangeable Applicators

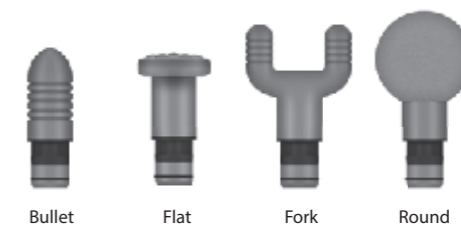USAGE INSTRUCTIONS

PLEASE READ PRIOR TO USE

CHARGING

- 1) Connect USB Type-C cable to the charging port.
- 2) The charging speed depends on what adapter is used. It is suggested to use a USB Type-C adapter supporting QC9V/18W, PD9V/18W, or PD15V/30W to fully charge in approximately 3-4.5 hours.
- 3) During charging, the power indicator flashes in green slowly. When the green light stops flashing, the device is fully charged.
- 4) The power indicator shows the battery level:
 - a) Flashing Red Light: 0-10%
 - b) Solid Red Light: 10-25%
 - c) Solid Blue Light: 25-75%
 - d) Solid Green Light: 75- 100%

DISCLAIMER

Please do not use the ELEEELS percussive massage device, or any massage device without first obtaining approval from your doctor if you have any of the following conditions: pregnancy, neuropathy, retinal damage, pacemaker user, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements, surgical pins or plates, or you have any concerns about your physical health.

WARRANTY

IF WITHIN ONE YEAR FROM THE DATE OF PURCHASE, THIS PRODUCT FAILS DUE TO A DEFECT IN THE MATERIAL OR WORKMANSHIP, ELEEELS WILL REPAIR OR REPLACE THE PRODUCT; OR FAILED COMPONENT(S), FREE OF CHARGE

THIS WARRANTY EXCLUDES:

- (a) Damage caused by accident, abuse, or mishandling during transport;
- (b) Damage arising from use of units from unauthorized repair;
- (c) Malfunction of units that are not used in accordance with ELEEELS care instructions;
- (d) Damage for which the cost of repair exceeds the cost of the product;
- (e) Deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises, and
- (f) Any product for which a dated proof of purchase cannot be provided.

Please visit below link for instruction manual in other languages:
<http://www.eleeeels.com/im>

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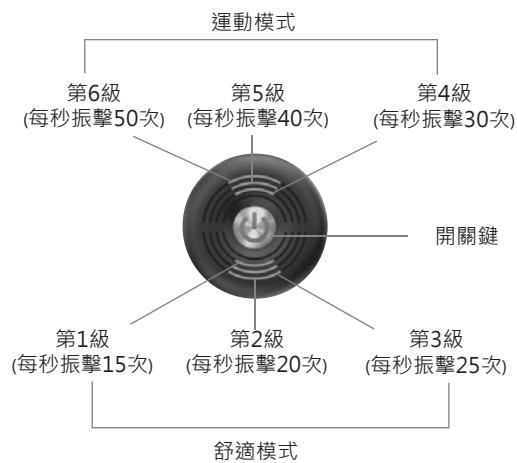
THIS WARRANTY EXCLUDES:

- (a) Damage caused by accident, abuse, or mishandling during transport;
- (b) Damage arising from use of units from unauthorized repair;
- (c) Malfunction of units that are not used in accordance with ELEEEELS care instructions;
- (d) Damage for which the cost of repair exceeds the cost of the product;
- (e) Deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises, and
- (f) Any product for which a dated proof of purchase cannot be provided.

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